


Product Information

APROTININ

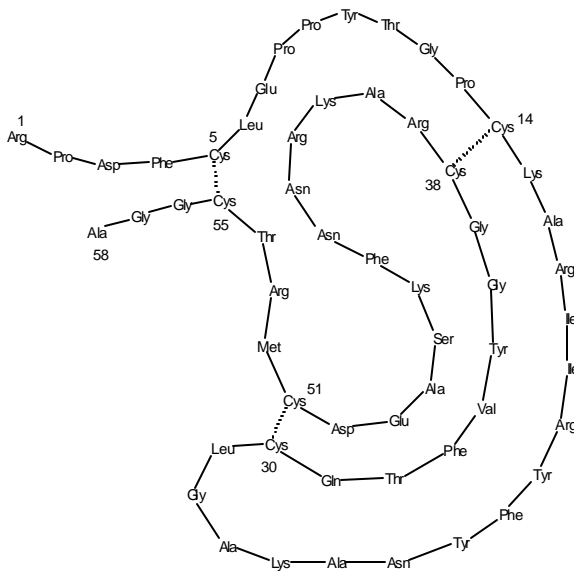

STABILITY / STORAGE AS SUPPLIED:

Products A1153, A4529, A3428, A6012 and A6279 should be stored at 2-8°C.

SOLUBILITY / SOLUTION STABILITY:

Aprotinin is freely soluble in water (>10 mg/mL) and in aqueous buffers of low ionic strengths.^{5,6} Dilute solutions are generally less stable than concentrated ones. Solution stability also depends on pH; values of 1-12 can be tolerated.⁴ Repeated freeze-thaw cycles should be avoided. The Cys¹⁴-Cys³⁸ disulfide bridge is readily split by reducing agents like β-mercaptoethanol.⁴ Due to its compact tertiary structure, aprotinin is relatively stable against denaturation due to high temperature, acids, alkalis, organic solvents or proteolytic degradation (only thermolysin has been found capable of degrading aprotinin after heating to 60-80°C).⁴ The high basicity of aprotinin causes it to adhere to commonly used dialysis tubing and even gel filtration matrices, but the use of acetylated materials and concentrated salt solutions (e.g., ≥0.1 M NaCl in buffer) minimizes the problem.⁴ Sterilization may be achieved by filtration through a 0.2 μm filter.⁵

SOLVENT	CONCENTRATION	STORAGE TEMP.	% LOSS/TIME
Sterile water with 0.9% NaCl and 0.9% benzyl alcohol, pH 5.7-6.2	10 mg/mL	0-5°C	<4.3%/year ⁵
2.5% Trichloroacetic acid	N/A	80°C	No loss ²
pH <12.6	N/A	N/A	No loss observed after 24 hrs. ⁷
pH >12	N/A	N/A	Irreversibly denatured ⁸
pH 7-8	0.065-1.95 µg/mL	4°C	About 1 week ⁶
pH 7-8	0.065-1.95 µg/mL	-20°C	>6 months ⁶

APROTININ

PRODUCT DESCRIPTION:

SIGMA PRODUCT NUMBER	DESCRIPTION
A-1153	Aprotinin from Bovine Lung Lyophilized powder Activity: 3-8 TIU/mg solid
A-4529	Aprotinin from Bovine Lung Lyophilized powder Activity: 3-7 TIU/mg solid Starting material: A-1153 Affinity purified to remove trace impurities.
A-3428	Aprotinin from Bovine Lung, Cell Culture Tested Lyophilized powder Activity: 3-7 TIU/mg solid Starting material: A-1153
A-6012	Aprotinin from Bovine Lung Aseptically filled solution in 0.9% NaCl and 0.9% benzyl alcohol. Activity: 5-10 TIU/mL solution
A-6279	Aprotinin from Bovine Lung Aseptically filled solution in 0.9% NaCl and 0.9% benzyl alcohol. Activity: 3-7 TIU/mg protein; 5-10 TIU/ml solution Similar to A-6012, but produced by Sigma.

USAGE:

Aprotinin is a competitive serine protease inhibitor which forms stable complexes with and blocks the active sites of enzymes. The binding is reversible, and most aprotinin-protease complexes dissociate at pH >10 or <3.²

APROTININ

USAGE: (continued)

ENZYME - SOURCE - CONDITION	INHIBITION (K_i = Dissociation Constant)
Acrosin	Weak inhibition ⁶
Chymotrypsin	$K_i = 9 \text{ nM}^9$
Chymotrypsinogen (bovine), pH 8.0	$K_i = 9 \text{ nM}^4$
CMP-N-Acetylneuraminate lactosylceramide α -2,3-sialyltransferase	74% Inhibition at 300 nm^9
Elastase (human leukocytes), pH 8.0	$K_i = 3.5 \text{ }\mu\text{M}^4$
Kallikrein (pancreatic), pH 8.0	$K_i = 1.0 \text{ nM}^4$
Kallikrein (plasma)	$K_i = 30 \text{ nM}; 100 \text{ nM}^9$
Kallikrein (tissue)	$K_i = 1 \text{ nM}^9$
Kallikrein (urine)	$K_i = 1.7 \text{ nM}^9$
Plasmin (porcine), pH 7.8	$K_i = 4.0 \text{ nM}^4$
Plasminogen activator	$K_i = 8 \text{ }\mu\text{M}; 27 \text{ }\mu\text{M}^9$
Trypsin (bovine), pH 8.0	$K_i = 0.06 \text{ pM}^4$
Trypsinogen (bovine), pH 8.0	$K_i = 1.8 \text{ }\mu\text{M}^4$
Tryptase TL-2	16% Inhibition at $10 \text{ }\mu\text{M}^9$
Urokinase (human), pH 8.8	$K_i = 8.0 \text{ }\mu\text{M}^4$

UNIT DEFINITION:

One Trypsin Inhibitor Unit (TIU) will decrease the activity of 2 trypsin units by 50%, where 1 trypsin unit will hydrolyze 1.0 μmole of N α -benzoyl-DL-arginine p-nitroanilide (BAPNA) per minute at pH 7.8 and 25°C. Another commonly used unit of activity is the KIU (Kallikrein Inhibitor Unit). A conversion factor for Aprotinin is: 1 TIU .1,300 KIU.⁵ A published ratio is: 1 TIU .1,025 KIU.¹⁰

APROTININ

REFERENCES:

1. *Merck Index*, 12th Ed., S. Budavari, Ed., # 796, p. 128 (1996).
2. *J. Gen. Physiol.*, 19, 991 (1936).
3. *Hoppe-Seyler's Z. Physiol. Chem.*, 192,1 (1930).
4. *Drug Res.*, 33(1), No. 4, 479 (1983).
5. Sigma data.
6. *Biochemica Information*, 1st Ed., J. Keesey, Ed., Boehringer Mannheim Biochemicals, p. 111, Indianapolis (1987).
7. *Biochem.*, 7, 3634 (1968).
8. *Life Sci.*, 28, 1861 (1981).
9. *Handbook of Enzyme Inhibitors*, 2nd Ed., Part B, H. Zollner, Ed., p. 572, VCH Verlagsgesellschaft, Weinheim (1993).
10. *Biotechnology*, p. 565 (June 1990).

Sigma warrants that its products conform to the information contained in this and other Sigma-Aldrich publications. Purchaser must determine the suitability of the product(s) for their particular use. Additional terms and conditions may apply. Please see reverse side of the invoice or packing slip.